



USER MANUAL

B60P

File: 2010-10-BP-131 B-50_0 GB

Table of contents:

1. Basic description	4
2. Completion	4
3. Technical data	5
4. Keys and indicators	6
5. Security principles	7
6. Environment protection	7
7. Preparation to work.....	8
8. Start-up	8
9. General rules	9
10. Balance checking	10
11. Scale adjustment.....	10
12. Connecting to a computer or a printer (optional)	10
13. Basic function.....	14
13.1 Normal weighing.....	14
13.2 Weighing with tare.....	15
13.3 Increased readability.....	15
14. Special functions description.....	16
14.1 Menu customization function (ACTIV and dEFAULT).....	17
14.2 Autozeroing function (AutoZE).....	18
14.3 Pieces counting (PCS)	19
14.4 Printer cooperation settings (PrInt).....	20
14.5 Serial port parameters setting (Port).....	21
14.6 Constant tare (tArE).....	22
14.7 Average calculation function (AVErAG)	24
14.8 Recipe weighing function (RECIPE).....	25
14.9 Percentage weighing function (PErC).....	26
14.10 Anti-disturbance filter function (FILtEr)	27
14.11 Function for maximum value indication (UP)	28
14.12 Force measuring function (nEWto)	29
14.13 Function for weighing large animals (LOC).....	30
14.14 Checkweighing function (thr)	31
14.15 Total weight function (totAL)	34

14.16	<i>Function for setting date and time (dAtE)</i>	36
14.17	<i>Function for selecting label number (LAbEL)</i>	37
14.18	<i>Setting time of stabilization function (Stb)</i>	38
14.19	<i>Entering reference zero function (Zero).....</i>	39
14.20	<i>Automatic switching off the scale function (AutoOF)</i>	40
14.21	<i>Density determining function (dEnSlT).....</i>	41
14.22	<i>Calculator for good packaging control (tP)</i>	43
14.23	<i>Statistical calculations function (StAt)</i>	46
14.24	<i>Paperweight calculation (PAP).....</i>	48
15.	<i>Maintenance and small repairs</i>	49
	<i>Declaration of Conformity</i>	50

1. Basic description

B60P platform scales are designed for general use. Thanks to built-in accumulator and a case to transport the scale (on demand) B60P is a great portable balance.

All scales are metrological tested - calibration or legal verification on demand.

Scales have following verification features:

- a seal protecting scale casing against opening,
- notified body stamps and green metrological marking placed on the balance name plate.

Renewing of scale legal verification is required when protective seal is violated or after period of 3 years starting from 1st December of year when first legal verification was performed.

In order to renew legal verification please contact authorized service of AXIS.

NACE classification: 29.24.23.

Certificates:



EC Type
Approval certificate
No. PL 04 022



Axis management System
Certificate DIN EN ISO 9001:2000
No. 78 100 6386

2. Completion


Basic set includes:

1. Scale
2. 4 legs
3. Feeder
4. User manual
5. Guarantee
6. Case to transport the scale (option on demand)

3. *Technical data*

Scale type	B60P
Maximal load	60kg
Minima load	400g
Readout graduation (d)	20g
Verification graduation (e)	20g
Accuracy class	III
Working temperature	-10 ÷ +40°C
Tare range	-60kg
Weighing time	<3s
Pan dimensions	400x500mm
Scale base dimensions	400x600mm
Scale height	80mm
Supply	AC: 220V, 50Hz DC: 12V, 1,5A
Internal supply	Lead-acid battery 6V/12Ah
Accumulator working time	about. 10 hours.
EC verification	✓

4. Keys and indicators

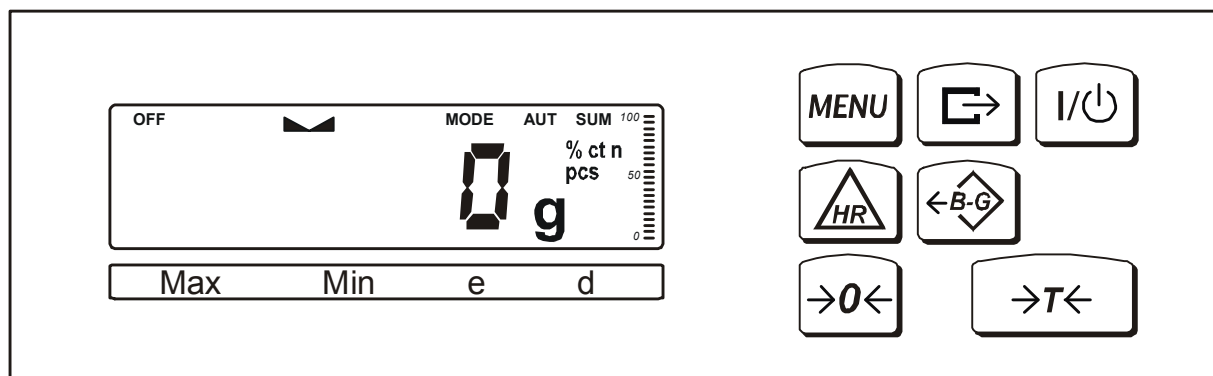


key	I/⏻	- switch-on / switch-off (standby),
key	→T←	- tare (subtract package weight from weighed mass)
key	B/G	- gross weigh indication switch,
key	→0←	- zero,
key	MENU	- special functions menu,
key	🖨️	- print-out,
key	HR	- increased indication resolution,
indicator	→0←	- zero indicator
indicator	▬	- result stabilisation indicator
indicator	NET	- net weight indicator (indication with subtracted tare)
indicator	MODE	- special function setting,
bar indicator		- total load indicator (graduated 0-100%)
indicator	OFF	- standby,
indicator	B/G	- gross mass (key B/G),
indicator	pcs	- pieces counting
Max, Min, d, e, III		- metrological parameters and accuracy class.


Attention:

→0←, B/G, HR keys and →0←, B/G, NET indicators don't work in scales without legalization with graduation $d \neq e$.

5. Security principles

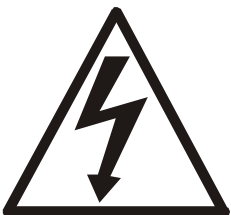	<p>To avoid electrical shock or damage of the scale or connected peripheral devices, it is necessary to follow the security rules below.</p>
<ul style="list-style-type: none"> • All repairs and necessary regulations can be made by authorised personnel only. • To avoid fire risk use a feeder of an appropriate type (supplied with the scale). Pay attention that supply voltage is compatible with specified technical data. • Do not use the scale when its cover is opened. • Do not use the scale in explosive conditions. • Do not use the scale in high humidity. • If the scale seems not to operate properly, unplug it from the mains and do not use until checked by authorised service. 	

6. Environment protection


	<p>According to legal regulations it is forbidden to dispose wasted electronic equipment in waste containers.</p>
<ul style="list-style-type: none"> • Please return wasted scale to the point of purchase or other company specialised in recycling of wasted electronic components. 	

7. Preparation to work

1. Take the scale out of the package removing protective foils.
2. Remove the lock nuts and unscrew lock screws (replace them with scale legs).
3. Place the scale on a stable ground not affected by mechanical vibrations and airflows.



Correct

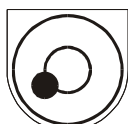

Wrong

4. Level the scale with the rotating rear legs so that the air bubble in the water-level at the back of the scale is in the middle.
5. Put on the pan.

6. Plug the feeder connector to the power socket at the back of the scale and plug the feeder to the mains. After self-tests and result stabilisation, the scale displays zero indication.
7. In scales with external feeder connect the feeder joint to supply socket situated on display column (sometimes situated on the cable). Connect the feeder to 230V.

8. Start-up

Leave the pan empty, plug the feeder to the mains (~230V/50Hz) and plug the feeder connector to the 12V power socket at the back of the scale.



The scale proceeds with following start-up actions:




Display test.



Basic electronic components tests: *C1*, *C-2*, ... *C-6*.

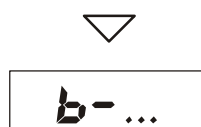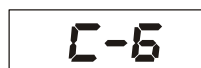

Program version.

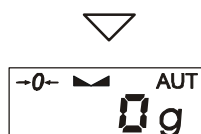

The scale is now ready to work.

9. General rules

1. Before each measurement make sure that zero indicator is displayed. If zero indicator does not displayed or “----“ communicate appears, press →0← key and wait until zero indication and zero indicator appears.
2. The scale is equipped with a tare equal to its range. To tare the scale press →T← key (left or right). Storing a tare value does not extend measuring range, but only subtracts it from a load placed on a pan. To make weight control easier and to avoid range overdrawing, the scale is equipped with a load indicator (graduated in percentages).
3. Weighing result should be read when the indicator "└ ┘" lights, which signalises stabilisation of a result.
4. When the scale is not used but it is necessary to be ready to work immediately, it can be switched off by pressing I/⏻ key. The scale reading system is then switched off to "standby" mode (signalled by the indicator "OFF" in version with the LCD display). To switch the scale on press I/⏻ key. The scale is immediately ready to operate maximum accuracy (after self tests).
5. Weighed sample should be placed in the centre of the pan.



Place the scale on a platform to avoid dropping weighed objects on the pan.



Do not overload the scale more then 20% of maximum load (Max).

6. Protect the scale against dust, aggressive dusts and liquids. To clean the scale wash it with water with soap and dry it afterwards.

10. Balance checking

It is advised to check scale indication accuracy before and after series of measurement using any load with known weight.


To check the scale with legal verification use a calibration weight with valid calibration certificate. In case permissible error is exceeded it is advised to contact the nearest service to calibrate the scale.

11. Scale adjustment



To adjust a balance it is necessary to break protective marks, therefore to calibrate the scale please contact the nearest service.

12. Connecting to a computer or a printer (optional)

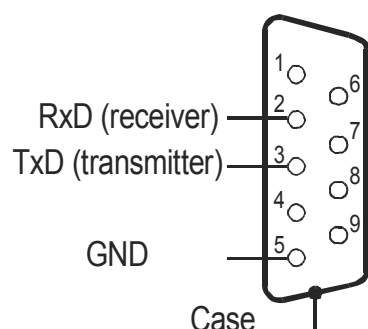



The scale is equipped with RS232C or RS485 (optional), which can be used to connect external devices such as a computer or a printer. Detailed way of communication data is defined by one of protocol:


- *LonG* – protocol to printer and computer
- *ELtron* – protocol to label printer

When cooperating with a printer data is send automatically after result stabilisation, but next transmission is possible after removing previously weighted sample.

RS485 (option)

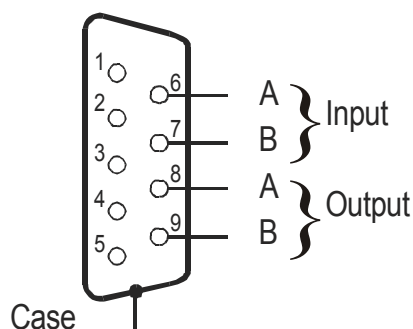

When protocol *Long* is chose scale sends result of weighting on an initialize signal from computer, or after pressed a  key. When automatic transmission mode is chose (function *LPt*), a printer data is send automatically after result stabilisation, but next transmission is possible after removing previously weighted sample. The scale sends following information: successive number of weighing, date and time (if clock module is installed)and weighing result. During transmission communicate Print is displayed

Scale with *Eltron* protocol sending record of data after pressed  key.

The scale sends following information: number of label, date and time (if clock module is installed) and weighing result. During transmission communicate Label is displayed.

In special features scale can be equipped in second RS232 interface for additional functions, ex. continues transmit of indicate to second display.

Computer must have a special program for cooperation with a data from a scale. This programs are procurement by Axis.

Detailed protocol description

LonG protocol

Communication parameters: 8 bits, 1 stop bit, no parity, baud rate 4800bps,

- initialising signal (data send order) – correspond to press  key

Computer→Scale: S I CR LF (53h 49h 0Dh 0Ah),

Scale→Computer: scale response according to description below (16 bytes):

Byte 1	- sign „-” or space
Byte 2	- space
Byte 3÷4	- digit or space
Byte 5÷9	- digit, decimal piont or space
Byte 10	- digit
Byte 11	- space
Byte 12	- k, l, c, p or space
Byte 13	- g, b, t, c or %
Byte 14	- space
Byte 15	- CR
Byte 16	- LF


Attention:

Network number different than zero (F..rS / *nr* function) changes scale working mode: communication with a computer is possible after logging the scale in with 02h scale number command. To log the scale out use 03h command.

For example: Using a program to test RS232 interface (program is available in [www.axis.pl / programy komputerowe](http://www.axis.pl/_programy_komputerowe)) for scale number 1 please write: \$0201 to log in, then *SI*, and write: \$03 to close communication.

- Scale tarring (calling $\rightarrow T \leftarrow$ key press) :
Computer \rightarrow Scale: S T CR LF (53h 54h 0Dh 0Ah),
Scale \rightarrow Computer: without response,
- Scale zeroing (calling $\rightarrow 0 \leftarrow$ key press):
Computer \rightarrow Scale: S Z CR LF (53h 5Ah 0Dh 0Ah),
Scale \rightarrow Computer: without response,
- Scale turning on / off (calling $\uparrow \Phi$ key press):
Computer \rightarrow Scale: S S CR LF (53h 53h 0Dh 0Ah),
Scale \rightarrow Computer: without response,
- Entering to special function menu (calling *MENU* key press):
Computer \rightarrow Scale: S F CR LF (53h 46h 0Dh 0Ah),
Scale \rightarrow Computer: without response,
- Setting low threshold value (option):
Computer \rightarrow Scale: S L *DI...DN* CR LF (53h 4Ch *DI...DN* 0Dh 0Ah)
DI...DN – threshold value, maximum 8 characters („-” – negative value, digits, dot – decimal separator), number of digits after dot should be the same as on scale display,
Scale \rightarrow Computer: without response,
- Example:
 - in order to set low threshold 1000g in scale B1.5 (d=0.5g) the following order should be sent:
S L 1 0 0 0 . 0 CR LF (53h 4Ch 31h 30h 30h 2Eh 30h 0Dh 0Ah),
 - in order to set low threshold 100kg in scale B150 (d=50g) the following order should be sent:
S L 1 0 0 . 0 0 CR LF (53h 4Ch 31h 30h 30h 2Eh 30h 30h 0Dh 0Ah),),
- Setting high threshold value (option):
Computer \rightarrow Scale: S H *DI...DN* CR LF (53h 48h *DI...DN* 0Dh 0Ah),
DI...DN – threshold value (see)
Scale \rightarrow Computer: without response.

Connecting cable WK-1 (scale – computer / 9-pin interface):

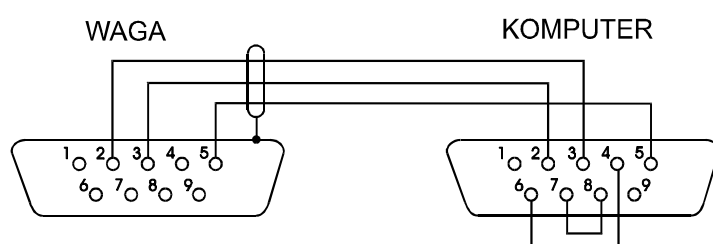

Protocol ELTRON description

Communication parameters: 8bits, 1stop, no parity, 9600bps,

- After use  key scale send:

Scale→Label printer : EPL-2 language instruction to initialize print of label:


US	- control instruction
FR"0001"	- number of labels
?	- starting instruction
hh:mm	- 5 signs: hours:minutes
yyyy.mm.dd	- 10 signs: year.month.day
mass + unit	- 10 signs: scale indicate + mass unit
P1	- ending instruction

Attention:


1. It is possible to place a constants sign (company name, product name).
2. The default number of label's is 0001 (label number 1).Different numbers of label are possible by using a *LABEL* function.
3. Label formula must be saved in printer – label forma must be designed in computer (by Zebra Designer program) and saved in printer.
Parameters and protocol of transmission must be compatible with printer label

13. Basic function

To make clear how to manage with each function, in further part of instruction descriptions are replaced with pictures.

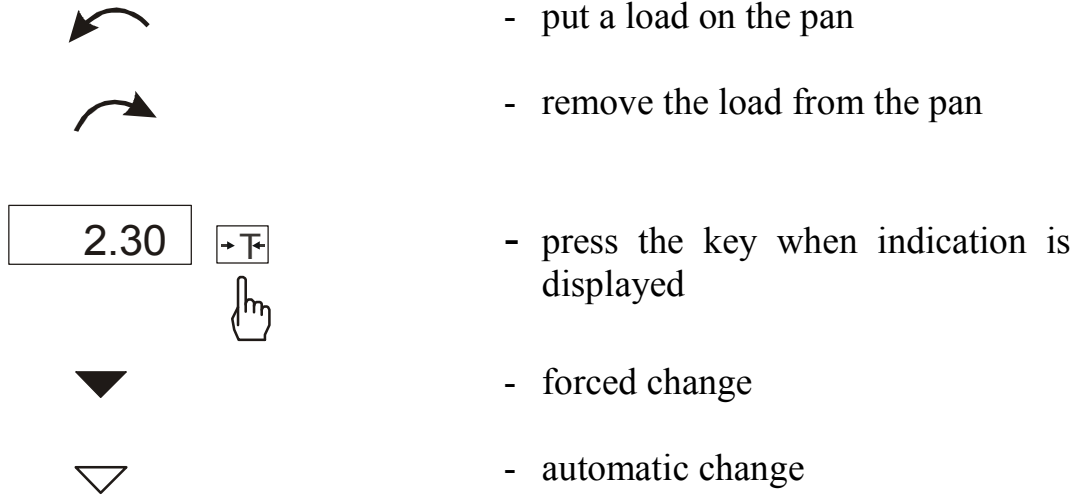


13.1 Normal weighing

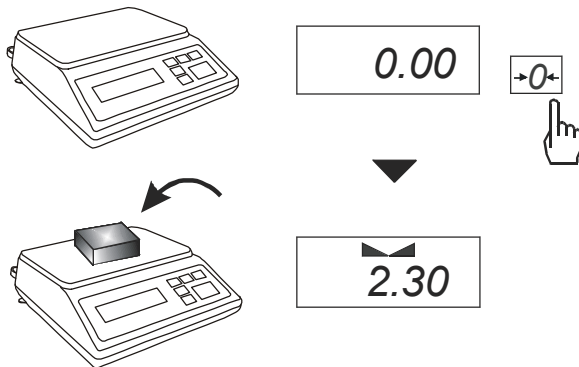

Press $\rightarrow 0 \leftarrow$ key ($\rightarrow T \leftarrow$ key in non-legalized scales), which zeros the scale, operates only when the pan is empty .

Weighing result should be read when the indicator " $\blacktriangle \blacktriangle$ " lights.

13.2 Weighing with tare

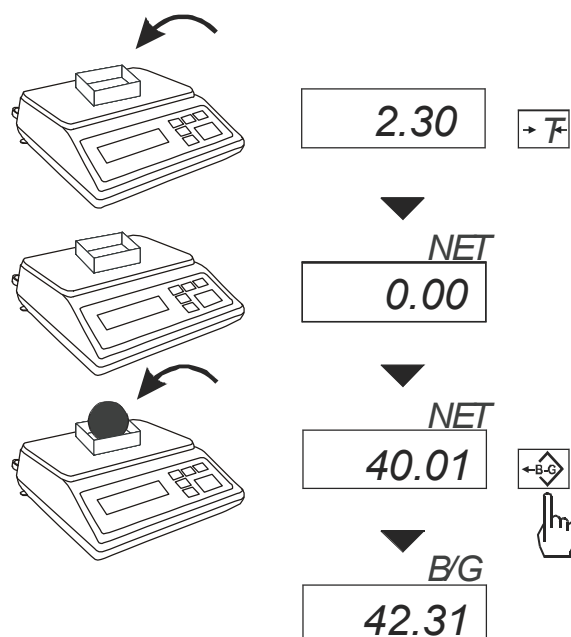


The scale is equipped with tare equal to its range.

Joint value tare and mass net can not cross a maximum of scale.

To display gross weight press *B/G* key.

Note:
Press *B/G* key to return to net weight indication.

13.3 Increased readability

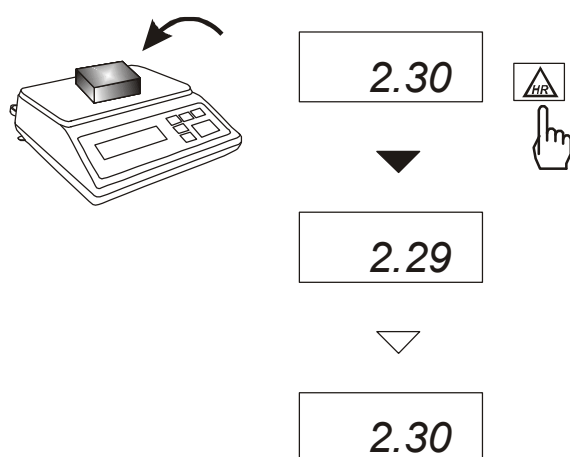

Press *HR* key to display the weighing result (for 5s.) with the highest readability possible. This function is especially helpful in scales with legal verification with $d=e$.

The weighing result with increased readability can be used for informational purposes only and cannot be printed or sent to a computer with **☞** key.

14. Special functions description

List of available functions:

- ❑ menu customization function (*ACtIV*),
- ❑ removal of all function from menu (*dEFAUL*).
- ❑ autozeroing function (*AutoZE*),
- ❑ pieces counting function (*PCS*),
- ❑ function for setting serial port working (*PrInt*)
- ❑ function for setting serial port (*Port*)
- ❑ entering tare function (*tArE*)
- ❑ recipe weighing function (*rECIPE*),
- ❑ weighing large animals function (*LOC*)
- ❑ force measuring function (*nEWto*)
- ❑ maximum value indication function (*UP*)
- ❑ anti-disturbance filter function (*FILtEr*)
- ❑ average calculating function (*AVErA*)
- ❑ percentage weighing function (*PErC*),
- ❑ extended calibration function (*CALib*)
- ❑ setting time of stabilisation function (*Stb*)
- ❑ selecting label number function (*LabEL*)
- ❑ automatic switching off scale function (*AutoOF*)
- ❑ entering reference zero (*Zero*)
- ❑ determining solids and liquids density function (*dEnSI*t)
- ❑ calculator for good packaging control (*tP*) - option*
- ❑ statistical calculations (*StAt*)- option*
- ❑ paperweight calculation function (*PAP*) - option*


function with additional equipment require:

- options with the clock:
 - setting current date and time function (*dAtE*)
 - total weight function (*totAL*)
- options with the transoptors connectors:
 - checkweighing function (*thr*)

* Functions offered with special version of scale software (with limited possibility of using other special functions).

User create own menu by choosing function in *ACtIV* function (described in chapter 14.1).

14.1 Menu customization function (*ACTIV* and *dEFAULT*)

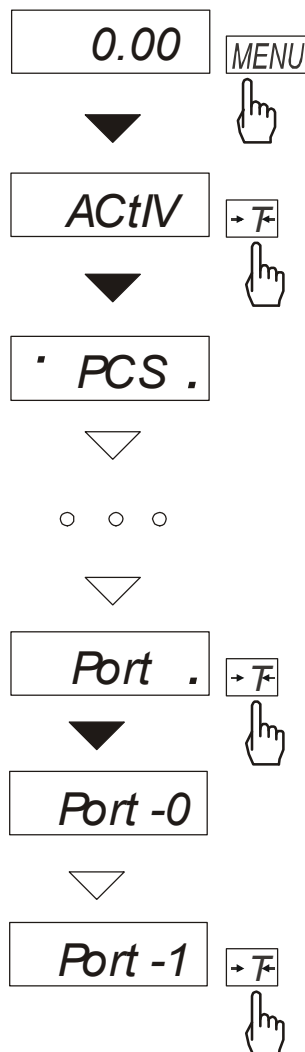

Among available user functions it is possible to select these, which should be displayed after pressing *MENU* key. It allows avoiding displaying whole list of available functions, which makes operation time longer.


Operation sequence shown in the pictures on the left causes adding function for setting serial interface RS232C parameters (*Port*) to function menu.

After switching on *ACTIV* function a dot is displayed on the right side (to distinguish from regular menu). Chosen functions are displayed with a dot on the left side.

In every moment, it is possible to restore primary (manufacture) settings choosing *dEFAULT* special function.

In order to remove function from menu in the last operation in place of selecting *Port -1* choose *Port -0*.

14.2 Autozeroing function (AutoZE)

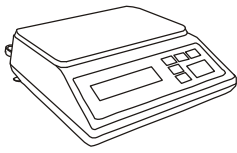

0.01 MENU



AutoZE →T←




Aut-0



Aut-1 →T←



▪ 0.00

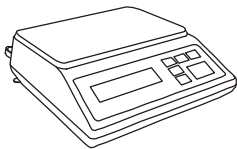


When *F..Aut* function is activated, the scale automatically ensures stable zero indication if the pan is empty or if zero indication was acquired by pressing *→T←* key.

To leave the function press *MENU* key, then with *→T←* key chose *AutoZE* and *Aut-0*.

Note:

Autozeroing function is activated automatically for 10 min. after switching-on.

14.3 Pieces counting (PCS)

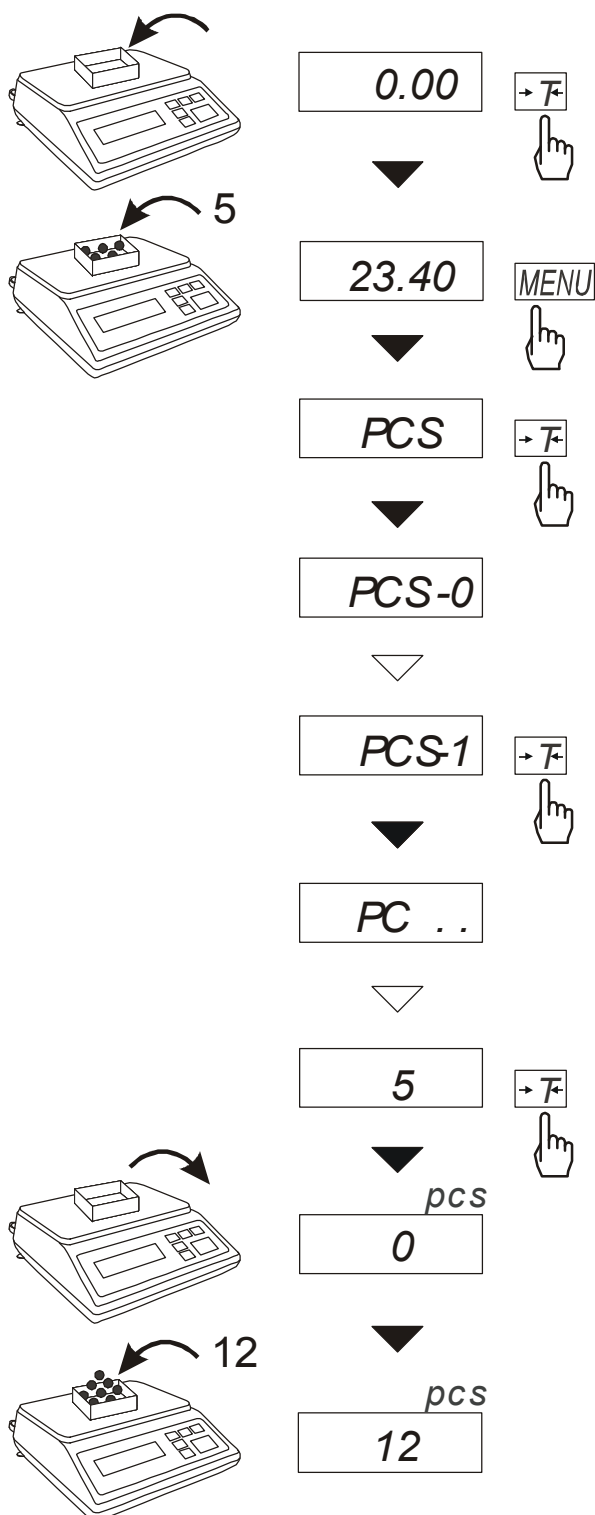

This function enables to count identical pieces, e.g. turnbuckles or buttons.

A measurement is performed in two phases:

- first phase - single piece weight calculation on the basis of defined pieces amount (5, 10, 20, 50, 100, 200 or 500 pieces),
- second phase – pieces counting.







It is advised that single piece weight is not less than one reading unit and sample weight used in first phase is bigger than 100 reading units.

To leave function press **MENU** key and then with **→T←** key chose **PCS** and **PCS-0**.

Note:

1. *Err-3* communicate signalises that a sample was not put on the pan. The same communicate appears if single piece weight is less than one reading unit (it is possible to count pieces but measuring error is bigger).
2. To chose previously used pieces amount select " _ _ " in first phase (in case no value was chosen, error communicate appears
3. During pieces counting **→T←** key function does not change.
4. In scales equipped with LCD display, weighing unit is visible and "□" sign is replaced with "pcs".

14.4 Printer cooperation settings (*PrInt*)

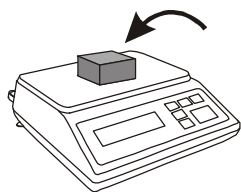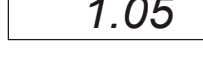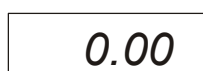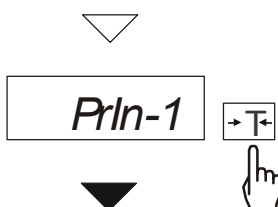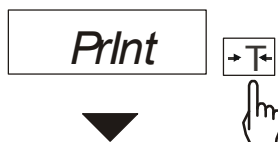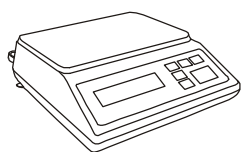


Activate the function for automatic serial port working mode (cooperation with a printer).


After activation the scale prints a header. Weighing result with a successive measurement number is printed automatically after result stabilisation (without using key).

To select computer cooperation mode (key activated and weighing results without successive numbers) press *MENU* key, then with $\rightarrow T \leftarrow$ key chose *PrIn-0* and *PrIn-1*.

14.5 Serial port parameters setting (Port)

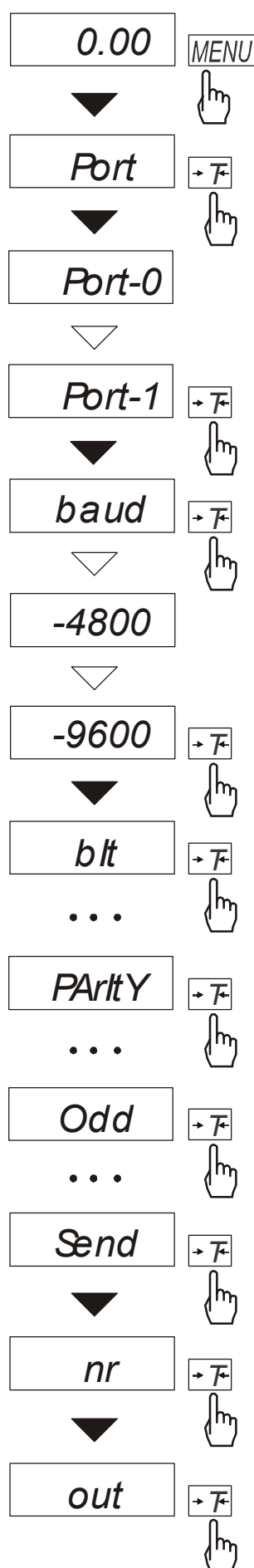

The function enables to set the following transmission parameters (standard parameters underlined:

- transmission protocol (*Prot*):
LonG - printer,
ELtron – label printer,
- transmission speed (*bAud*: 1200, 4800, 9600, ...),
- the number of bits in a byte (*bit*: 7, 8),
- parity control (*PARtY*: 0, 1; *Odd*: 0, 1),
- network number when working in multistand computer system (when working as a single scale the value should be “0”),
- continuous transmission – without using  key, approx. 10 results per second (*SEnd*: 0, 1).

Protocol *Eltron* automatically activated function *LABEL*.


To set desired transmission parameters activate *Port* function, choose appropriate parameter and press $\rightarrow T \leftarrow$ key to accept needed parameter value. The example at the left presents how to set transmission speed value to 9600bps.

To leave the function choose *out* option.

14.6 Constant tare (tArE)

This function enables to measure gross weight of a sample placed in a container of a known weigh value (stored in the memory) and to display calculated net weight of the sample. Tare value is recalled from the memory with $\rightarrow 0 \leftarrow$ key when the pan is empty. Tare value may be entered using the keypad or by sampling container weight from the pan.

Operation sequence:

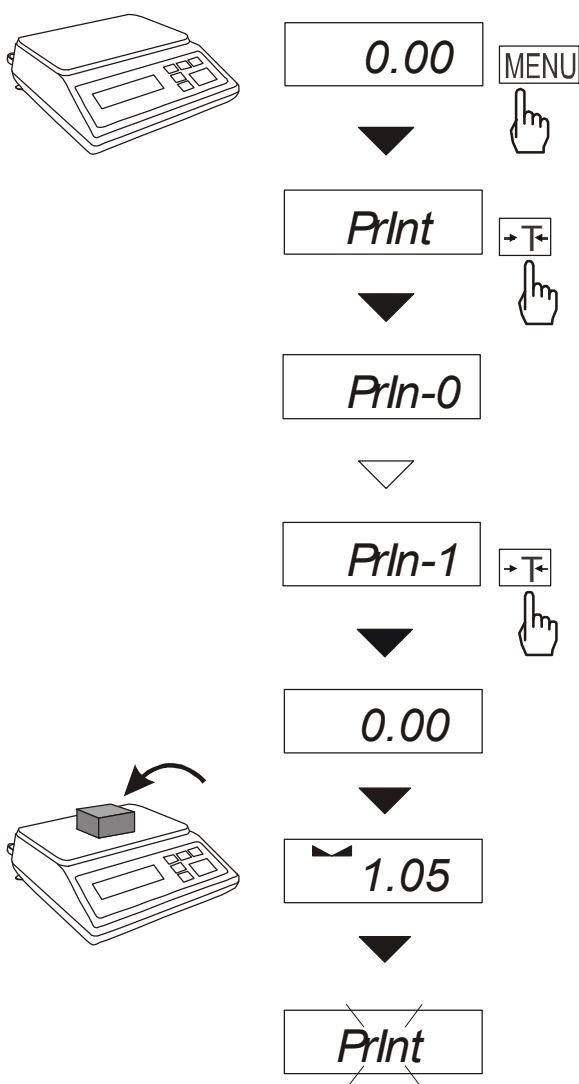

The following options are possible:

- *tAr-0* – leave the function,
- *tAr-1* – activate the function with the previous tare value,
- *tAr-2* – sample tare value from the pan,
- *tAr-3* – enter tare value with keys: $\rightarrow 0 \leftarrow$, $\rightarrow T \leftarrow$ and MENU
- *tAr-4* – printout a setting value of tare


If the function is active, *NET* indicator is displayed.

Options *tAr-1* enables to activate the function with previous tare value after leaving the function with *tAr-0* option.

Note:

Tare value is stored in memory also after unplugging the scale from the mains.

Weighing with constant tare:

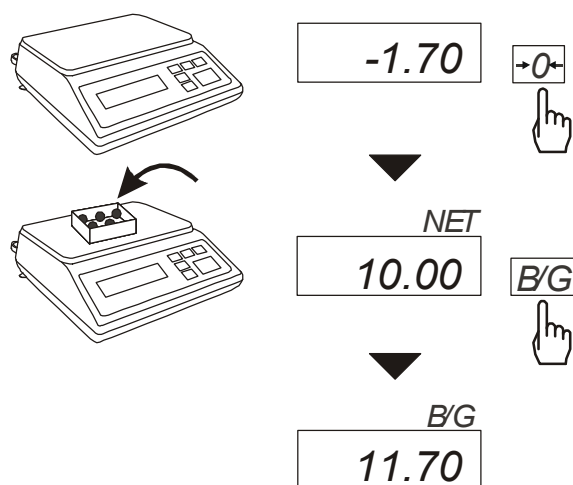

When *tAr* function is activated, press **→0←** key to zero the indication and to recall tare value from the memory. Tare value is displayed with ”-” sign.


B/G key enables for instant switching between net and gross weight.

Note:

When the pan is empty **→T←** key does not operate – to tare the scale use **→0←** key.

14.7 Average calculation function (AVERAG)

The function allows for calculating average value of performed measurement series. During series of measurements successive results are registered automatically when scale indication is stabilised.

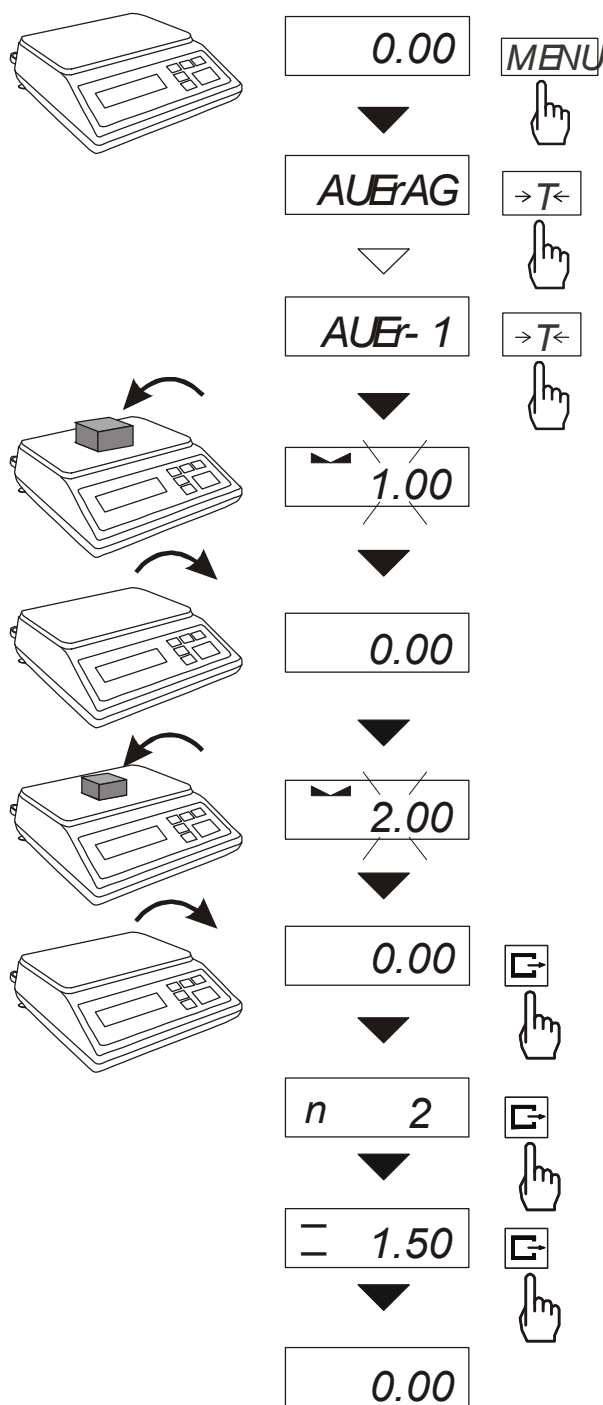

Press *MENU* key and select *AUErAG* pressing $\rightarrow T \leftarrow$ key.

The following will be shown successively on display:

- *AUEr-0* – leaving function,
- *AUEr-1* – measurement with average calculation..

Select *AUEr-1*. It will allow weighing with simultaneous storing results into summing register for average calculation

Measurement registration is performed automatically in the moment when scale indication becomes stable. Short time displaying of „- -” denotes that load can be taken off and new one can be put on. Results above scale Min are registered only. Number of measurements is limited to 9999.

In order to read average value \square key should be used.


- first pressing causes displaying number of measurements (n).
- second displaying causes displaying average value (=).
- third pressing allows continuing average calculation.

If a printer is connected to scale the following report will be printed :

Date: ...	Time. ...
MEASUREMENTS No	= ...
AVERAGE VALUE	= ...

In order to finish calculation press *MENU* key, and then select *AUEr* i *AUEr-0*.

14.8 Recipe weighing function (RECIPE)

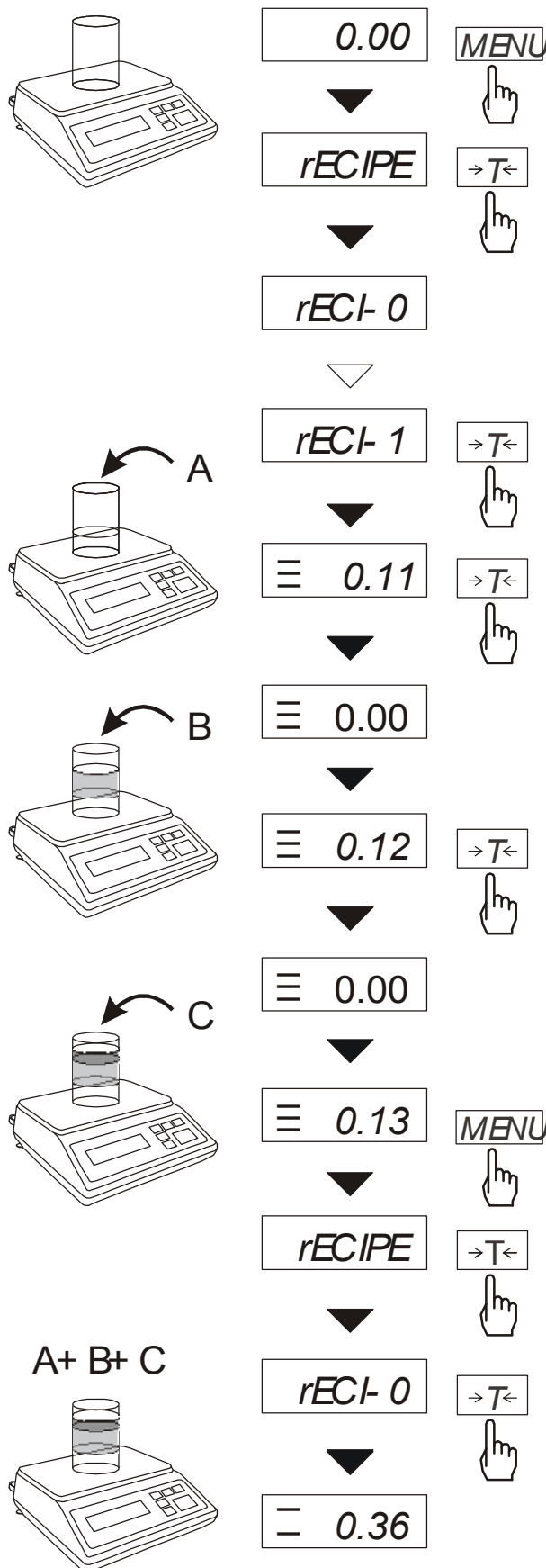

This function allows for separate weighing of several ingredients in one container with a possibility of control total weight of all weighed components.

The function has the following options:

- *rECI-0* – leave the function with possibility of reading to read total weight,
- *rECI-1* – start recipe weighing
- *rECI-2* – continue previous recipe.

When preparing a recipe successive ingredients (A, B, C, etc.) are weighed each time starting from zero indication. In order to allow this after weighing of each ingredient tare the scale.

After weighing of several ingredients reading total weight is possible (despite scale taring). In order to do that press *MENU* key, select *rECIPE* function once more and use *rECI-0* option.


Sign „≡” signals total weight indication. Recipe is finished when →T← key is pressed.

When „≡” sign is displayed recipe continuing is possible. *rECI-2* option is used for that.

Note:

Sign „≡” on the left side of display informs about *rECIPE* function activity.

14.9 Percentage weighing function (PErC)

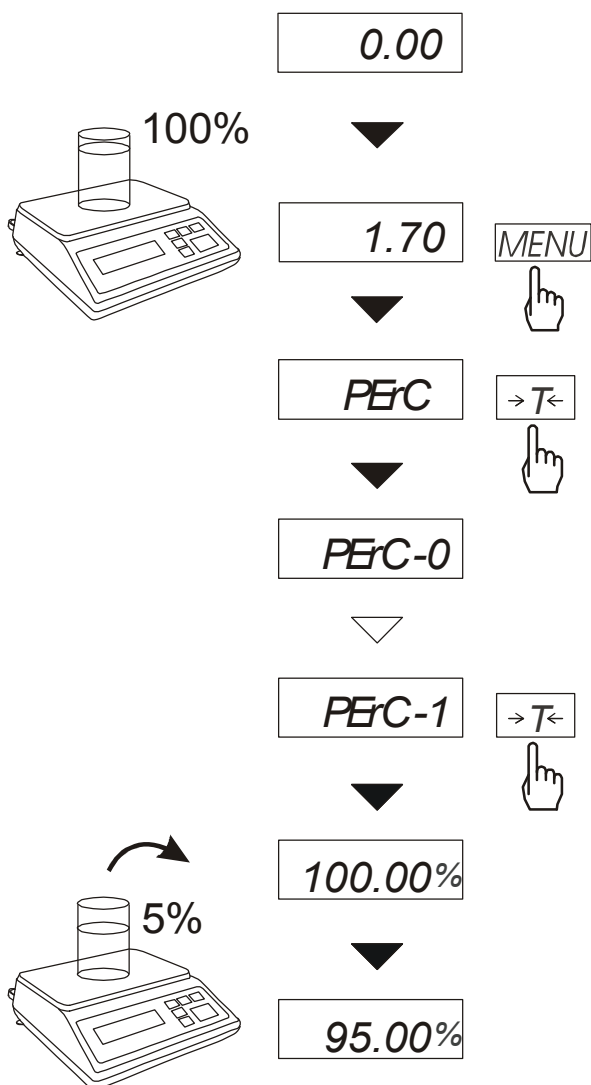

This function allows displaying weighing result in percents.

A measurement is performed in two phases:

- first phase – weighing a reference sample (100%),
- second phase – measuring specific sample as a percentage of the reference sample.

Weighing result is displayed in different format, depending on the weight value of reference sample. For weight values of reference sample $0 \div 3,5\%$ of weighing range result is displayed in format *100*, for range $3,5 \div 35\%$ - in format *100.0*, and above 35% - in format *100.00*.

„%” sign is replaced with „ $\frac{\square}{\square}$ ” indicator.


The function has the following options:

- *PErC-0* – disable the function,
- *PErC-1* – set current scale indication as 100% and activate percentage weighing,
- *PErC-2* – continue percentage weighing after exiting to normal weighing.

Note:

During percentage weighing **→T←** key has its normal function.

14.10 Anti-disturbance filter function (*FILtEr*)

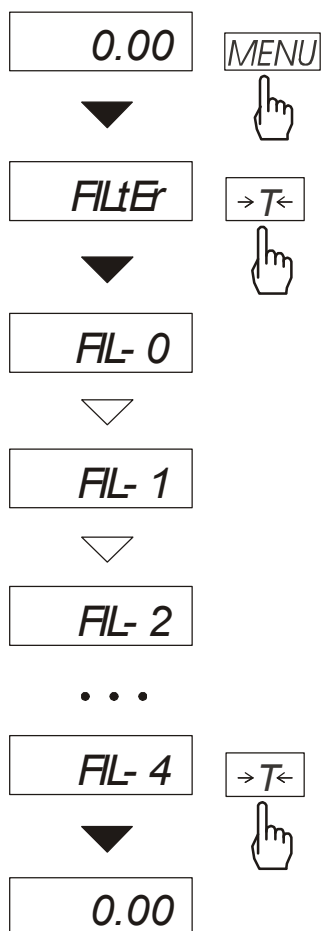

This function allows using digital filter with selected intensity during weighing. Filter reduces the influence of mechanical vibrations (air blasts, base vibrations) on measurement result.

Press *MENU* key and select *FILtEr* pressing **→T←** key.

The following options will be shown successively on display:


- *FIL-0* – work without
- *FIL-1* - filter I (weak)
- *FIL-2* - filter II (medium)
- *FIL-3* - filter III (sharp)
- *FIL-4* - filter IV (very sharp)

Select one of four filters. This will cause starting weighing with selected filter.

In order to go back to normal weighing use *MENU* key once more and choose *FILtEr* and *FIL-0*.

14.11 Function for maximum value indication (UP)

This function allows holding on display maximum value shown by the scale in a while.

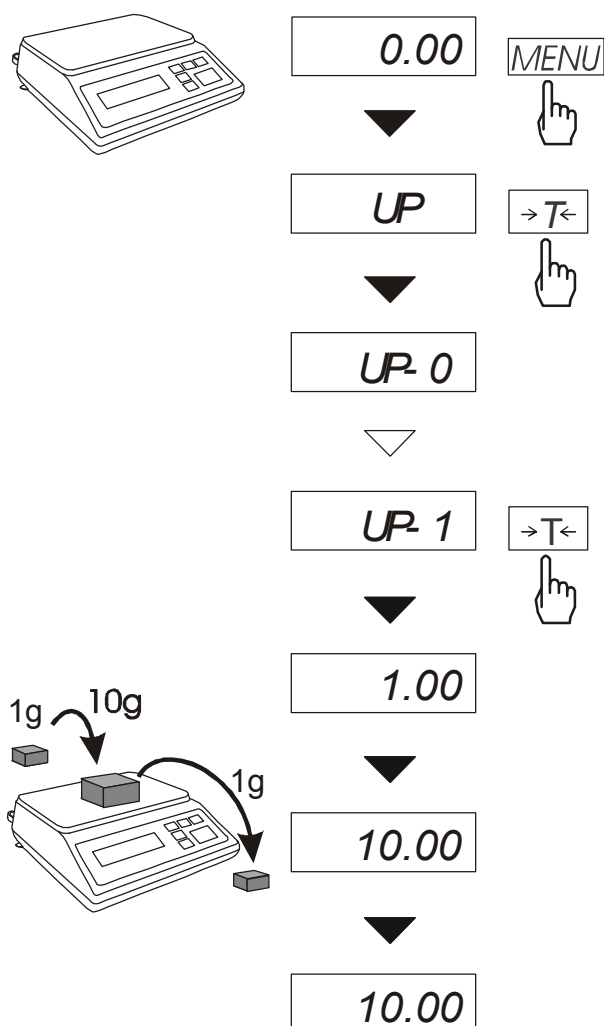

Before measurement scale should be tared.


After using *MENU* key and selecting *UP* function the highest mass result will be hold on display.

Pressing $\rightarrow T \leftarrow$ key will cause result zeroing.

Note:

Autozeroing function and the stabilisation indicator are deactivated when UP function is running. Weighing result is continuously averaged from 5 measurements.

14.12 Force measuring function (*nEWto*)

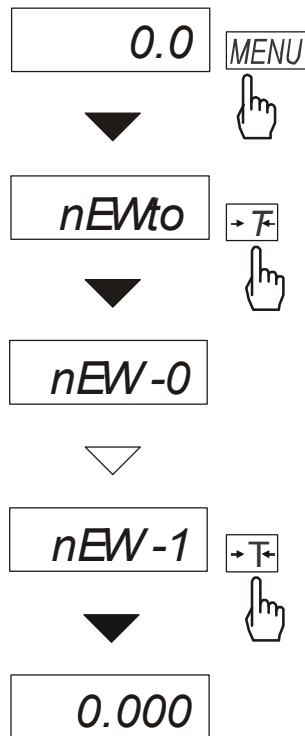

Function activation causes displaying result in force units (N).



Press *MENU* key.

Using **→T←** key choose *NEWto* function, and then *NEW-1*.

Note: 1N≈0,1019kg

14.13 Function for weighing large animals (LOC)

The function allows weighing animal moving on the scale.

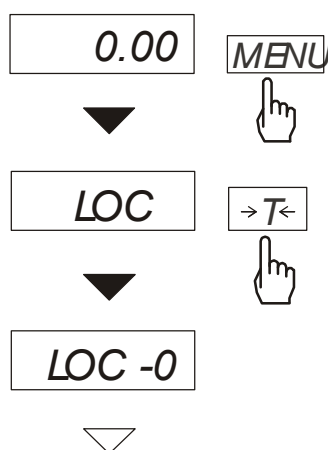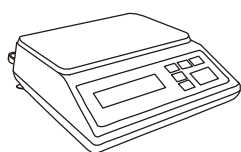

Press *MENU* key.

When *LOC* function is displayed press $\rightarrow T \leftarrow$ key.



The following options appear on display successively:

- *LOC-0* – leave the function,
- *LOC-1* – automatic weighing after loading the scale,
- *LOC-2* – the measurement initiated manually by pressing $\rightarrow T \leftarrow$ key.



When *LOC-1F..-1* is displayed press $\rightarrow T \leftarrow$ key.

Tare the scale using $\rightarrow T \leftarrow$ key if necessary and place the animal on the pan.

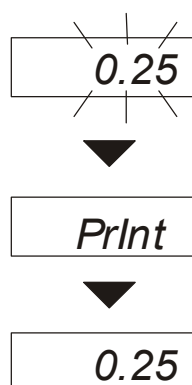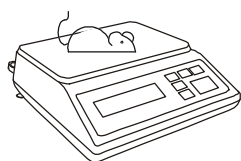

Wait until the weighing result is averaged – scale display will be blinking. Then scale will show stable averaged result and will send it through serial port. Final result is displayed on the display and send via serial port to computer or printer.

The result remains on display for about 30 second.

Important notes:

1. The loads less than Min are not averaged.
2. In the case when placing the animal takes more than 5s, it is advised to use *LOC-2* option (measurement initiated manually). It will allow performing measurement in right moment pressing $\rightarrow T \leftarrow$ key.

14.14 Checkweighing function (*thr*)

This function allows comparing weighing result with two programmed reference values: lower and upper threshold. Comparison result is signalled with indicators (MIN, OK, MAX) and sound signal generated when threshold values are exceeded. If comparison result is:

- smaller than lower threshold – the scale signals MIN (yellow colour),
- between threshold values - the scale signals OK (green colour, with the short sound signal),
- greater than upper threshold - the scale signals MAX (red colour, long sound signal).


The checkweighing results can be use to control:

- optical indicator (*Indication* mode),
- batching devices (*Batching* mode).

Standard scale is set for cooperation with optical indicator.


On outputs P1-P3 (*Relays* socket) short-circuit states appear as result of comparison scale indication with threshold values.

On the chart below output states are shown during increasing load on the scale for both working modes:

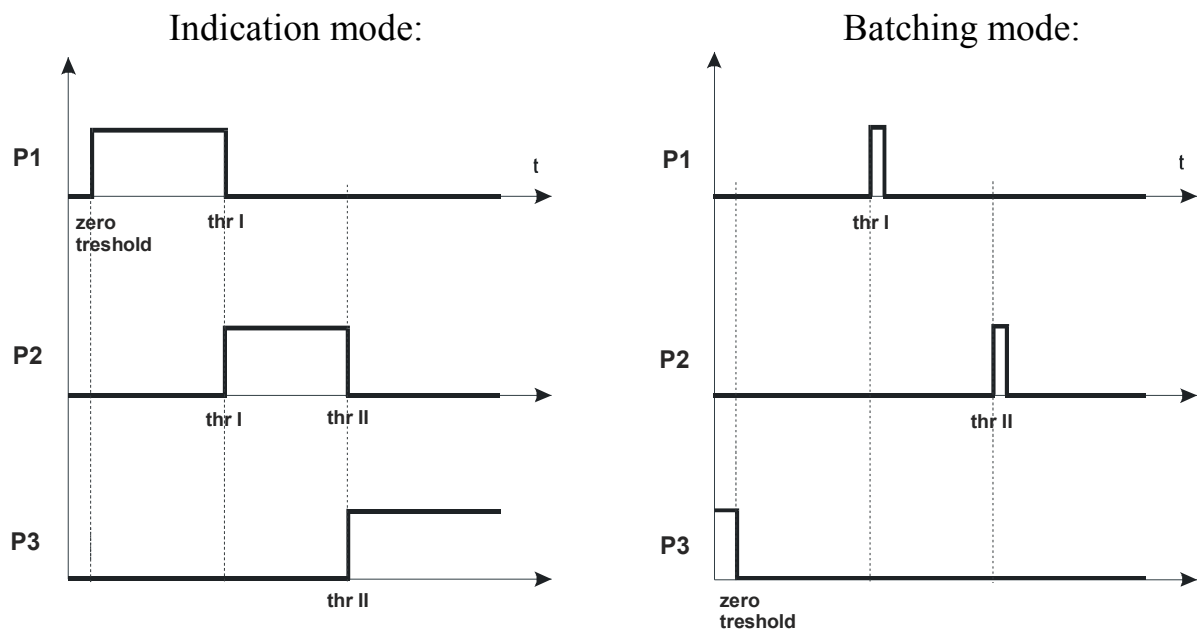

In *Batching* mode on P1 (thr I) and P2 (thr II) outputs short-circuit impulses appears for time of 0,5s. On P3 (zero) output short-circuit state appears when indication does not exceed threshold value signalling zero load.

Relays connection diagram:

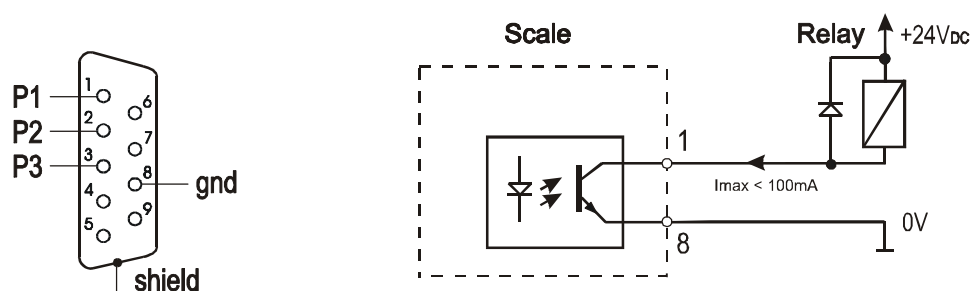




Relays output is the open collector transistor output with load capacity 100mA / 24V. Transmitter inputs must be protected with diodes, e.g. 1N4148.

It is advised to use MS3K/P electronic board (sold separately), consisting of RM96P transmitters, with DC24V input voltage and AC250V, 3A output.

Important notes:

1. After switching the scale on, both thresholds are set to maximum values.
2. When setting upper threshold value, pay attention that its value is not below lower threshold value.
3. Setting lower and upper threshold value is possible after sending appropriate orders from computer, what is described in scale user manual.

Operation sequence:

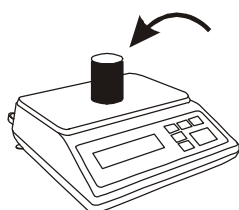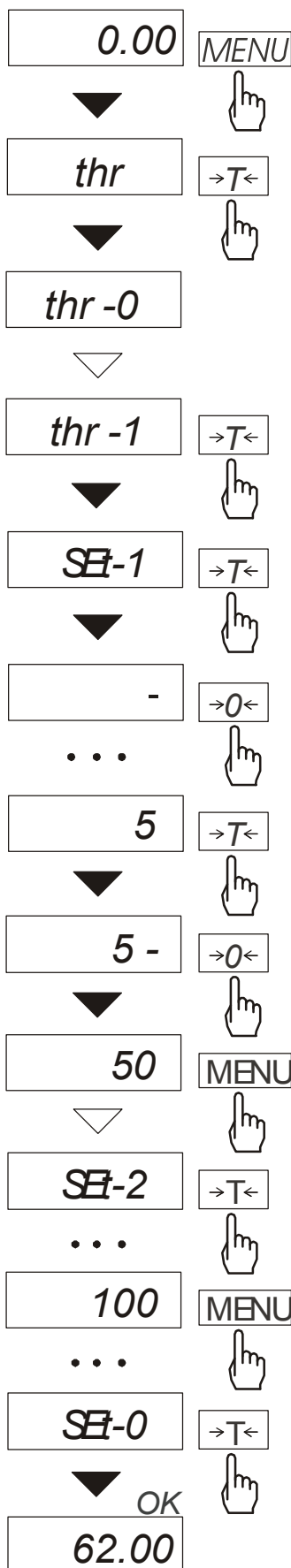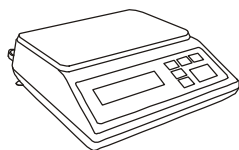

Press *MENU* key and choose *thr* pressing $\rightarrow T \leftarrow$ key.

The following options are displayed successively:

- *thr-0* – deactivate the function,
- *thr-1* – activate the function,
- *thr-P* – check last threshold values (press $\square \rightarrow$ key several times),
- *thr-t* – choose *Relays* socket mode:
 - 0 – exit to weighing
 - 1 – *Batching* mode
 - 2 – *Indication* mode.

Choose *thr-1* option using $\rightarrow T \leftarrow$ key. The following options for entering thresholds are displayed:

- *SEt-0* - go to weighing with signalling threshold excess,
- *SEt-1* - set lower threshold value,
- *SEt-2* - set upper threshold value,
- *SEt-3* - set zero signalisation threshold.

Using $\rightarrow T \leftarrow$ key select *SEt-1* option.

Set lower threshold value using the following keys:

- $\rightarrow 0 \leftarrow$ - digit increase,
- $\square \rightarrow$ - decimal point,
- $\rightarrow T \leftarrow$ - move to next digit,
- MENU* - finish.

Then select *SEt-2* option and enter upper threshold value.


Choosing *Set-0* option will cause starting work with signalisation of exceeding thresholds and zero.

To change *Relays* socket mode use *thr-t* option. Default option is *Indication*.

To leave the function, press *MENU* key and then choose *thr* and *thr-0* options.

14.15 Total weight function (totAL)

The function allows calculating total weight for series of measurements, which can be greater than scale capacity. It allows calculating total weight as well as average value.

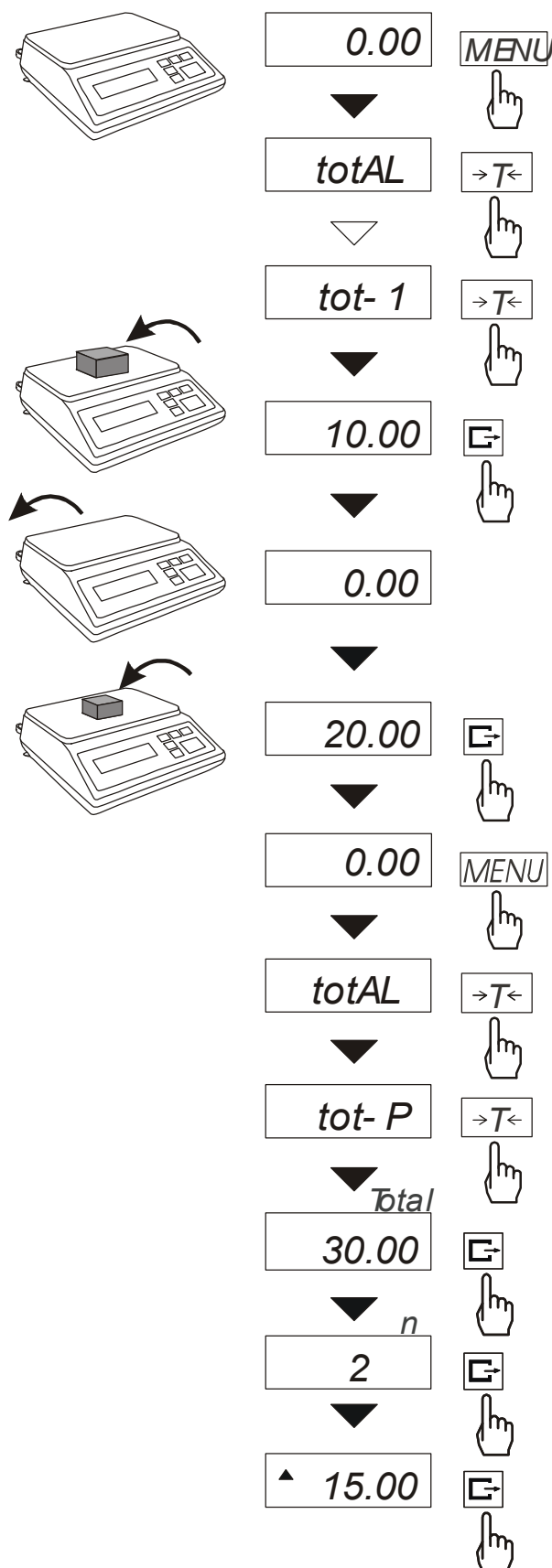

Press **MENU** key.

When **totAL** is displayed press **→T←** key.

The following options appear successively:

- **tot-P** - report printout without clearing total register,
- **tot-0** - clearing total register, report printout and leaving the function,
- **tot-1** - working with receipt printout after each measurement,
- **tot-2** - working without receipt printout.

Press **→T←** key when **tot-1** is displayed. Perform measurement series pressing **□→** key for storing results into total register.

In order to print and display results enter to function choosing **total** and **tot-P** option from menu

The results are display in the following sequence:

- total weight (\equiv)
 - number of registered measurements (n),
 - average value (\equiv),
- regarding that moving to display successive result is performed after pressing **□→** key.

In order to go back to total weighing without zeroing total register press **□→** key third time.

To leave the function with clearing total register, select *total* function from menu and choose *tot-0* option. When It will cause the scale prints the communicate informing about clearing registers.

The form of receipt after each measurement:

Date: ...	Time. ...
measurement no	weight
measurement no	weight

Report form:

Date: ...	Time. ...
TOTAL WEIGHT	=
NUMBER OF SAMPLES	=
AVERAGE VALUE	=

Note:

When the scale has not an internal clock, Date and Time do not appear on printout.

Maximum number of measurements 99 999.


Maximum total load 99 999 000d.

The weighing unit of the total value from the register (Total) is the same as the weighing unit stated on the keypad or is 1000 times greater, what is signalled by "o" indicator at the left of the display.

If the registered value is too big to be displayed, "E" communicate appears on the display.

If the number of series is too high and cannot be displayed, "Err1" communicate appears on the display.

14.16 Function for setting date and time (dAtE)

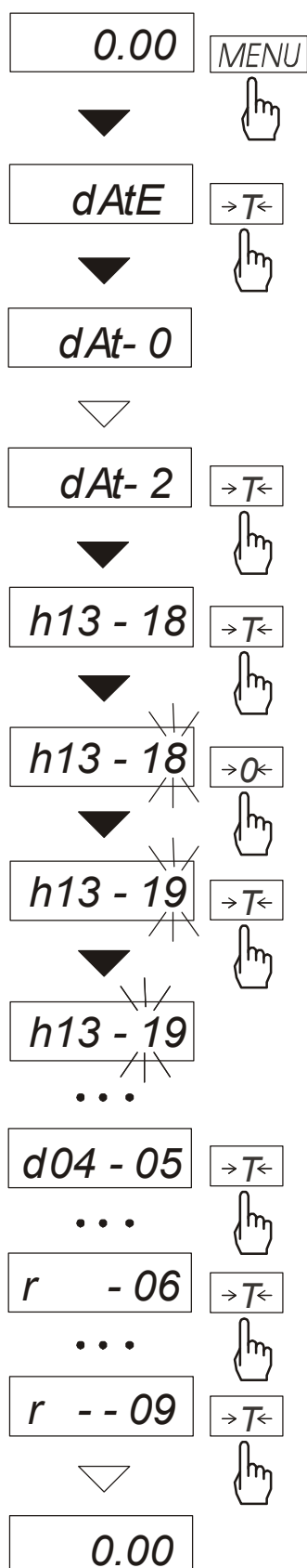

The function allows setting current date and time of scale internal clock and mode of its use.

The function has the following options:

- *dAt-0* – deactivate date and time during printout of current weighing result,
- *dAt-1* – activate date and time during printout of current indication (↵ key),
- *dAt-2* - change current date and time.

The example at the left presents how to set current date and time using *dAt-2* option.

On successive positions digits are changing automatically or manually using →0← key several times.

In order to choose appropriate digit and move to the next position use →T← key.

After setting proper date and time it should be activated with *dAt-1* option.


Time format: *h gg – mm*
(gg - hour, m - minute).

Date format: *d mm – dd*
(m - month, d - day).

Year format: *r - rr*
(r - two last year digits).

14.17 Function for selecting label number (LAbEL)

The function appears in scales with *ELTRON* data transmission protocol. This protocol allows printing scale indication and optionally date and time on label printer, as variable texts. Other data, e.g. company address, product name, its bar code can appear on label as constant fields. Label forms used by user, named as a numeric value (max. 4 digits) should be previously stored in printer memory according to printer user manual. Choosing label form is performed by entering label number using *LAbEL* function.

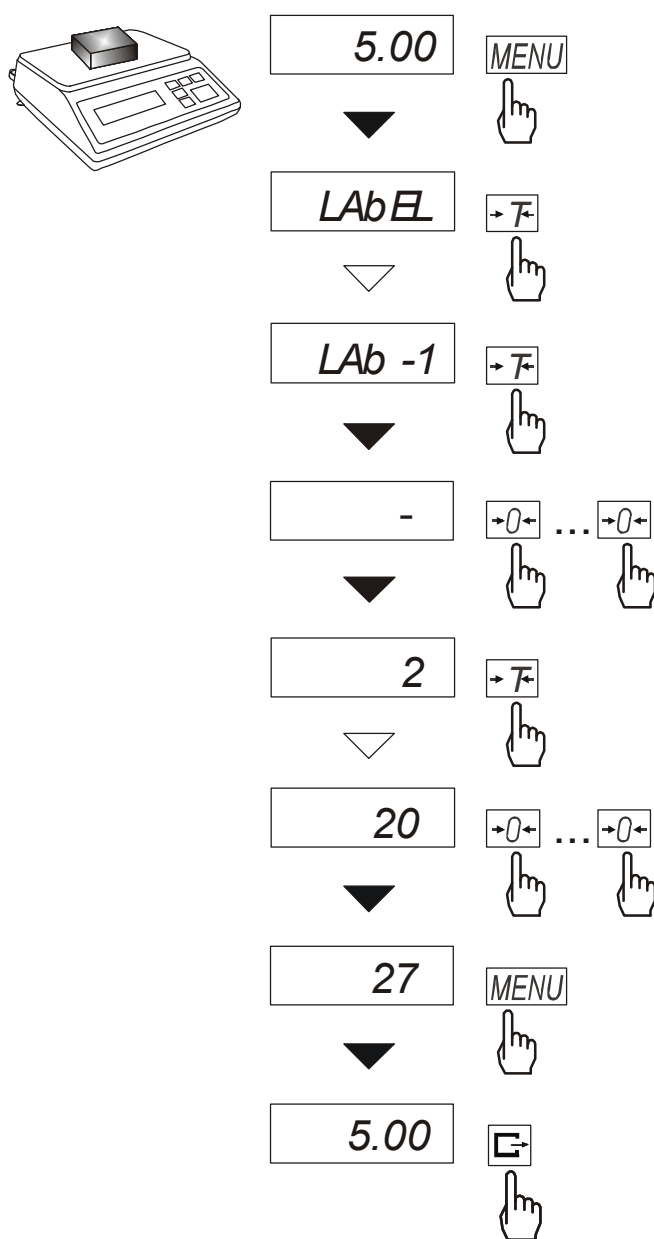

Press *MENU* key.

When *LAbEL* is displayed press **→T←** key.

The following options appear successively on display:

- *LAb-0* – leave without changes,
- *LAb-1* – enter label number,

Using **→T←** key select *LAb- 1*.

For entering label number the following keys should be used:

- 0←** – increase digit,
- T←** – move to next digit,
- MENU** – finish entering.


After putting load on and pressing **→T←** key data is sent to a printer.

Format of data sent to label printer:

US	(55 53 0D 0A)
FR"0001"	(46 52 22 30 30 30 31 22 0D 0A)
?	(3F 0D 0A)
00:00	(30 30 3A 30 30 0D 0A)
2000.00.00	(32 30 30 30 2E 30 30 2E 30 30 0D 0A)
10 g	(20 20 20 20 20 31 30 20 20 67 0D 0A)
P1	(50 31 0D 0A)

14.18 Setting time of stabilization function (Stb)

The function allows changing stabilisation time of scale indication and connected with it the time of waiting for starting result printout on a printer connected to the scale.

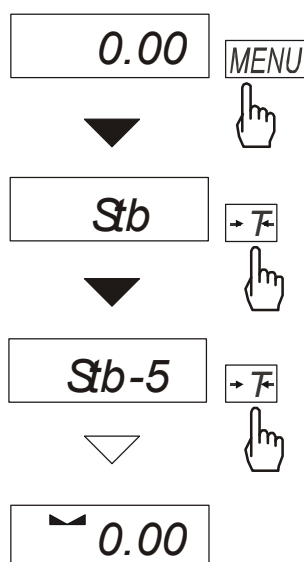

Press *MENU* key.

When *Stb* is displayed press **→T←** key.

The following options appear successively on display:

- *Stb-0* – deactivate the function,
- *Stb-1* – the longest stabilisation time,
- *Stb-2* – long stabilisation time,
- *Stb-3* – medium stabilisation time,
- *Stb-4* – shorter stabilisation time,
- *Stb-5* – the shortest stabilisation time.


After choosing one of the options weighing with activated filter is started..

In order to go back to normal work of the scale choose *Stb-0* option.

14.19 Entering reference zero function (Zero)

Note: This function is enabled in non legalized scales only.

ZER function allows entering new value of reference zero (value referred to empty pan) without need of contacting with authorised service centre.

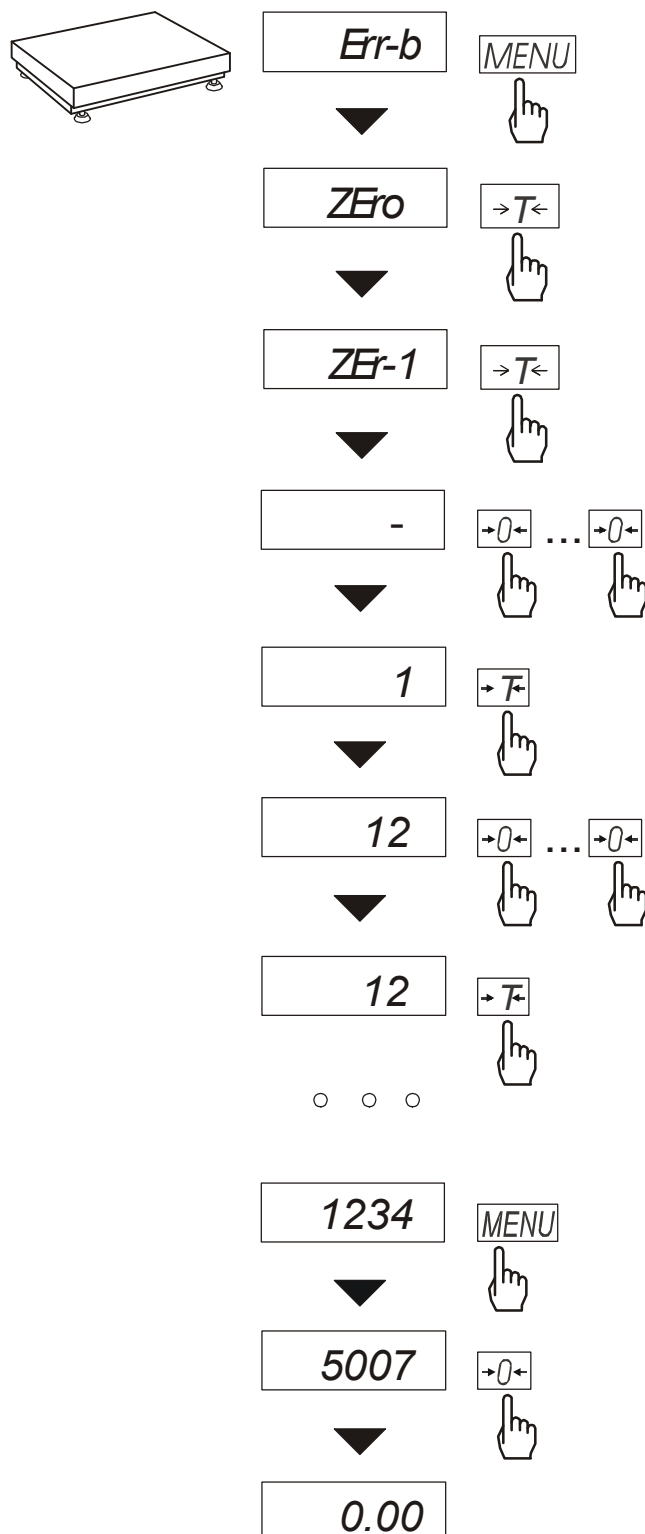

Press **MENU** key.

When **Zero** is displayed press **→T←** key.

The following options appear successively on display:

ZEr-0 – activate function,

ZEr-1 – enter new zero value,

ZEr-2 – enter new protecting code.

Using **→T←** key, choose **ZEr-1** and enter access code for function (in new scale it is 1234). Direct result from A/C converter will appear on scale display.


When the pan is empty press **→0←** key.

Wait for finishing zeroing process.

In order to change access code use **ZEr-2** option. Entering value is performed similarly as with **ZEr-1** option.

14.20 Automatic switching off the scale function (AutoOF)

The function is helpful in scales supplied from accumulator. The function causes scale to switch off automatically after c.a. 10 minutes of not using it. Switching function on causes last entered zero and tare values are remembered in scale memory. After next start-up of the scale these values are restored.

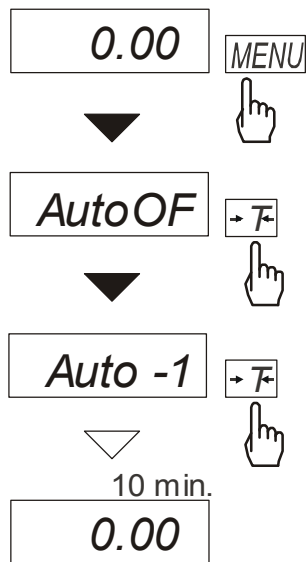

Press *MENU* key.

When *AutoOF* is displayed press **→T←** key.

The following options appear successively on display:


Auto-0 – activate function,

Auto-1 – deactivate function.

After function activation scale will control changes of its indications. If the scale will not be used and scale indications will not be changing, after time of c.a. 10 minutes the scale will switch off remembering its settings (zero and tare).

Switching on the scale is performed after pressing **I/⊙** key.

14.21 Density determining function (dEnSlT)

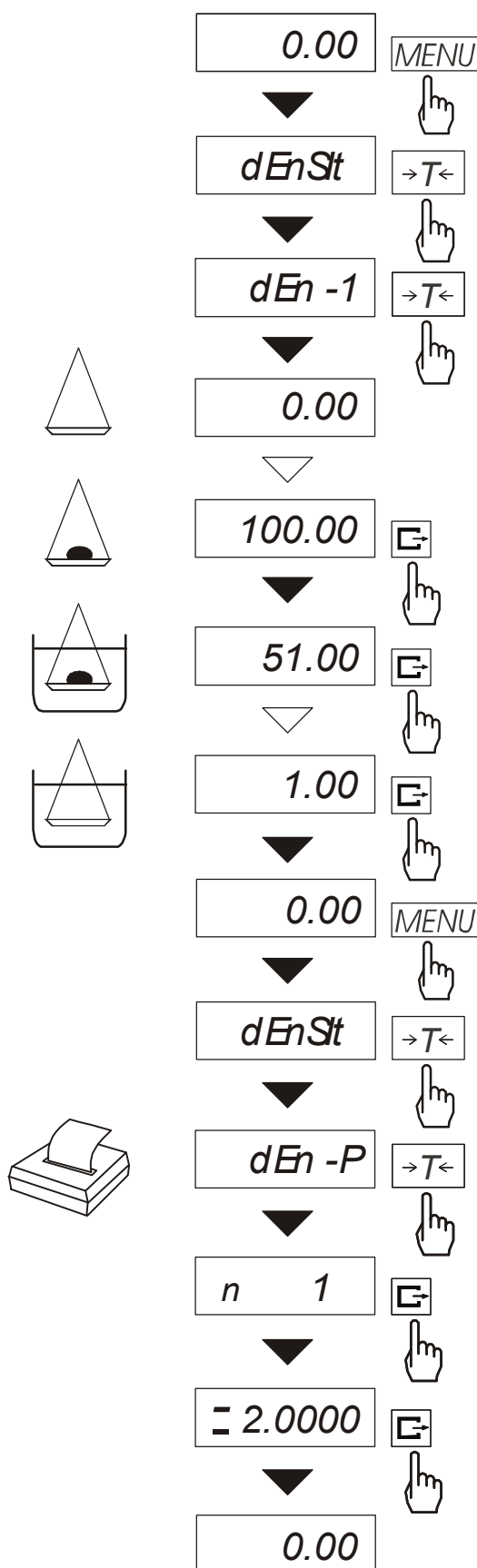

The function allows determining solid density basing on its weight in the air and in liquid of known density using the formula below:

$$g = \frac{m_1}{m_1 - m_2 + m_3} * g_{\text{liquid}}$$

where: m_1 -weight in the air

m_2 - weight in liquid

m_3 - hanger weight

g_{liquid} – liquid density

By default:

$$g_{\text{liquid}} = 1 \text{ g/cm}^3$$

(for distilled water).

When using liquid other than distilled water, choose *dEnSlT* from menu and use *dEn-2* option to enter liquid density taking into consideration its temperature.

To enter value use the following keys:

- 0←** - digit increase,
- - decimal point,
- T←** - move to next digit,
- MENU** - finish.

The measurement is performed in three phases:


- measurement in the air,
- measurement in liquid,
- hanger weighing

To read density determination result, enter *dEnSlT* function menu and use *dEn-P* option. First pressing **□** key causes displaying successive measurement number. Second pressing **□** key causes displaying and printing result, and then going to the next density measurement.


If a printer is connected to the scale, printout of solid density measurement results will be performed in the following form:

Date: ...	Time. ...
MEASUREMENT No.	=
WEIGHT in air	= g
WEIGHT in liquid	= g
HANGER WEIGHT	= g
LOAD DENSITY	= g/cm ³
LIQUID DENSITY	= g/cm ³

It is recommended to use pan hanging below balance, weighing in the air and in liquid is performed then in the following way.

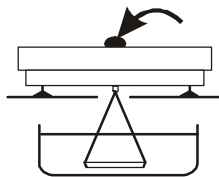

Phase I: measurement in the air.

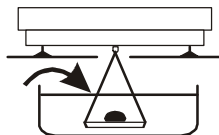

Phase II: measurement in liquid.

14.22 Calculator for good packaging control (tP)

Note: This function is enable for special orders substitutable with other functions.

The function allows performing not destroying control of 60, 100 or 160 samples of packed goods. Calculation formula complies with the description of reference method described in act of law from 26.07.2001 (with later changes).

Scale operator enters nominal net mass (Qn) and sample quantity taken for control (quantity refers to total quantity of samples 1 and 2 in description of reference method).

Accepting successive measured packages (storing in register) is performed automatically after load is put on and scale indication is stable.

Each time after load is put on printout is performed: measurement number, result, date and hour.

Accepting next measurement is possible after taking last load off.

For obtained measurement series scale calculates:

- \bar{x} -average mass as (sum x)/n

- Min -minimal mass in n samples

- Max -maximal mass in n samples

- R = Max-Min -the difference between max and min

- S -standard deviation
$$S = \sqrt{\frac{1}{(n-1)} \sum_n (x_n - \bar{x})^2}$$

- defective samples count -count of samples with mass $Qn-2T1 < x < Qn-T1$

-disqualifying samples count -count of samples with mass $< Qn-2T1$

Results of statistical calculations and histogram are printed on a printer.

Operation sequence:

1. Press *MENU* key.
 2. When "F..tP" is displayed press $\rightarrow T \leftarrow$ key.
- The following options appear successively on display:

- "F..0"- leave function,
- "F..1"- measurements,
- "F..2"- enter product data: product code and batch size,
- "F..3"- enter control parameters: nominal mass Qn and measurement number,
- "F..4"- select mode for data transmission (automatic or manual).


[Entering product data]

3. Press $\rightarrow T \leftarrow$ key when "F..-2" is displayed.
The following options appear successively on display:
 - "Cod"- product code [6 digits],
 - "n" - batch size,
 - "out" - exit,
4. Press $\rightarrow T \leftarrow$ key when "Cod" is displayed.
Previously entered product code will be displayed one by one with "out" option.
In order to enter new product key press $\rightarrow T \leftarrow$ key when previous code is displayed, and after "-" appears use keys:
 - $\rightarrow 0 \leftarrow$ - change digit,
 - $\rightarrow T \leftarrow$ - move to next digit,
 - MENU* - finish.
 If previous code is valid, choose "out" pressing $\rightarrow T \leftarrow$ key.
5. Press $\rightarrow T \leftarrow$ key when "n" is displayed.
Previously entered batch size will be displayed one by one with "out" option.
In order to enter new batch size press $\rightarrow T \leftarrow$ key when previous value is displayed using keys:
 - $\rightarrow 0 \leftarrow$, *T* and *MENU* as above.
6. Press $\rightarrow T \leftarrow$ key when "out" is displayed.

[Entering control parameters]

7. Press $\rightarrow T \leftarrow$ key when "F..-3" is displayed.
The following options appear successively on display:
 - "O_n"- nominal mass (Qn),
 - "-20", "-60", "-100", "-160", "n" – select available measurement quantity (total quantity of samples),
 - "out" - exit,
8. Press $\rightarrow T \leftarrow$ key when "O_n" is displayed.
Previously entered nominal mass will be displayed one by one with "out" option.
In order to enter new nominal mass press $\rightarrow T \leftarrow$ key when previous value is displayed using keys: $\rightarrow 0 \leftarrow$, *T* and *MENU* as above.
9. Select measurement quantity pressing $\rightarrow T \leftarrow$ key. Selected quantity should comply with requirements of good packaging act of law (it depends on control kind and batch size). Letter n denotes full control.
10. Press $\rightarrow T \leftarrow$ key when "out" is displayed.

[Selecting mode for data transmission]

11. Press $\rightarrow T \leftarrow$ key when "F..-4" is displayed.
The following options appear successively on display:
 - „Auto” – automatic data transmission after scale indication is stable,
 - „Recz” – manual data transmission by operator using  key,
 - „out” – exit.
12. Choose appropriate option pressing $\rightarrow T \leftarrow$ key.
13. Press $\rightarrow T \leftarrow$ key when "out" is displayed.

[Measurements]

14. Press *MENU* key.
15. When "F..-tP" is displayed press $\rightarrow T \leftarrow$ key. Select "F..-1".
16. Printout is performed.
17. Put successive good packages

Successive results will be printed in table with indication of their values using "*" character referred to limit values.

18. After performing last measurement "END" text will appear and summarising report of control results will be printed:

PLACE OF
CONTROL:

Date: ... Time: ...

CONTROL NO.:

SCALE TYPE : ...

FACTORY NO. : ...

BATCH SIZE : ...

VALUE Qn : ...

VALUE Qn-T1 : ...

VALUE Qn-2T1 : ...

	Qn-2T1	Qn	Qn+2T1
... g	*		
... g			*
.....			

19. In order to finish working with the function and reset results register, press *MENU* key, and when "F..tP" and "F..0" is displayed press $\rightarrow T \leftarrow$ key.

Note: Activating TP function causes that indicators signal exceeding limit values Qn-2T1 and Qn+2T1.

Date: ... Time: ...

MEASUREMENT COUNT = ...

QUALIFYING AVERAGE = ...

STANDARD DEVIAT. S = ...

HISTOGRAM

<Qn-2T1 -n_{2T1}

Qn-2T1 -----

A	n _A	█
B	n _B	██
C	n _C	███
D	n _D	██
E	n _E	████
F	n _F	█████
G	n _G	██████
H	n _H	█████
I	n _I	███
J	n _J	███
Qn-2T1		-----
>Qn-2T1 +n _{2T1}		█

RESULT : ...

CONTROLLER : _____

14.23 Statistical calculations function (StAt)

Attention: Function is available on demand and it replaces other special functions.

This function evaluates from series of measurements (max 500) statistical parameters of weighting process. Adding successively measurements to register is automatic and it occur after the scale is loaded and its indications stabilize.

After each loading printout is made with: number of measurements, result, date and time.

Next measurement is made after taking off earlier load.

For the obtained measurements series the scale evaluates:

- n - number of samples
- sum x - sum of all samples $sum_x = \sum x_n$
- \bar{x} -average value (sum x)/n
- min -minimal value from n samples
- max -maximal value from n samples
- max-min -maximal value minus minima value
- S -standard deviation $S = \sqrt{\frac{1}{(n-1)} \sum_n (x_n - \bar{x})^2}$
- srel -variance factor $srel = \frac{S}{\bar{x}}$

Statistical calculations results can be printed.

Order of operations:

1. Press *MENU* key.
2. When *F..StA* is displayed press $\rightarrow T \leftarrow$ key.
The following options are displayed:
 - *F..P* – statistical data printout,
 - *F..0* - out of function, register zeroing, statistic data printout,
 - *F..I* - enter or continue function,
3. Press $\rightarrow T \leftarrow$ key when *F..I* is displayed.
4. Put on successively objects on pan, (remove after indication stabilization) in order to add them to measurement register.
5. In order to obtain printed statistical results from measurements series press *MENU* key.
When sign *F..StA* is displayed, press $\rightarrow T \leftarrow$ key.
When *F..P* is displayed press again $\rightarrow T \leftarrow$ key.

This will cause printout of calculated statistics and histogram :

LSL - allowable lower value,

USL - allowable upper value,

A, B, C, .. - measurement intervals,

n_A ... - amount of measurements in A interval;
measurement is in A interval if it is bigger or equal to A interval threshold and smaller than B interval threshold.

n_B ... - amount of measurements in B interval;
measurement is in B interval if it is bigger or equal to B interval threshold and smaller than C interval threshold.

Thresholds are printed under histogram.

-NG - amount of measurements under allowable lower value

+NG - amount of measurements above allowable upper value


To finish work with this function and zeroing result register press *F* key, then during "*F..-StA*" and "*F..-0*" is displayed, press $\rightarrow T \leftarrow$ key. This will cause printing message about register zeroing.

Data:	...	Hour.	...
SAMPLES	=		...
TOTAL MASS	=		...
AVER MASS	=		...
MIN MASS	=		...
MAX MASS	=		...
MAX - MIN	=		...
S	=		...
SREL	=		...

HISTOGRAM

LSL ...
USL ...
DIV ...

-NG ...
LSL

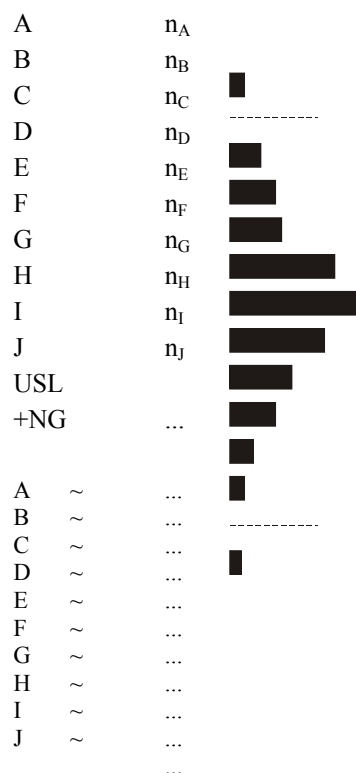

Measurement made by : _____


Cooperation between statistics function with computer and printer.

The scale can be equipped with two RS232C connections marked as RS232C-I (computer) and RS232C-II (printer). After each printer data printout, identical set of data is send to computer. After S A CR LF (53h 49h 0Dh 0Ah) initialization signal is sent by computer, the scale sends to computer statistical data contained in histogram.

14.24 Paperweight calculation (PAP)

Note: This function is enable for special orders substitutable with other functions.

This function enables to calculate paperweight of 1m^2 of paper basing on samples of known area. For quick access, the function is accessible directly by pressing *MENU* key.

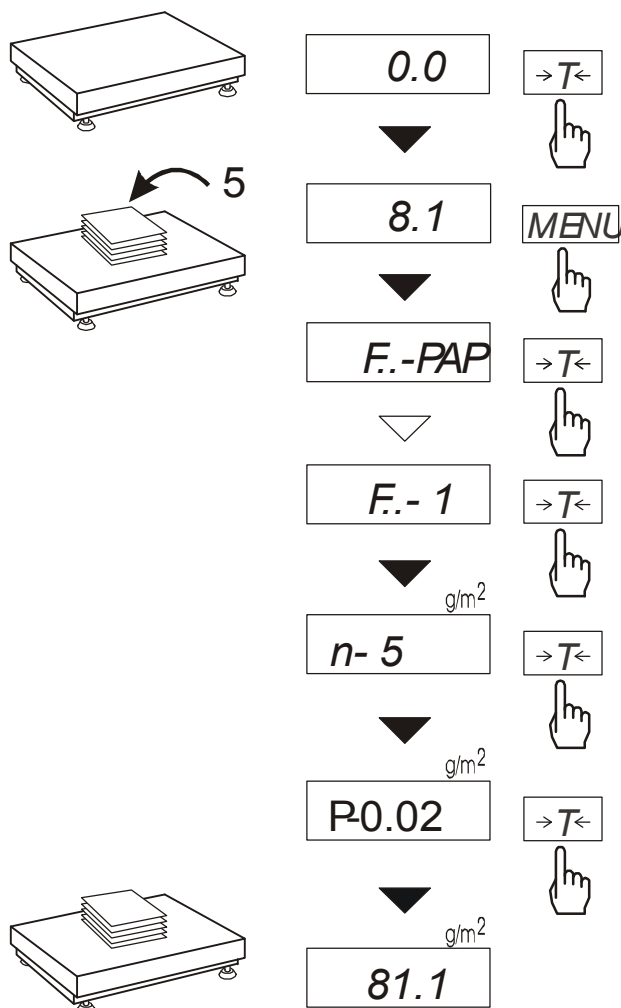

The balance must be tared just before the measurement.

Place the specific sample quantity of the same area (possible values: 1, 2, 5, 10, 20, 50, 100).

Press *MENU* key to access Function Menu. To enter the function press $\rightarrow T \leftarrow$ key when *F..-PAP* is displayed.

n is the number of samples placed on the pan. To choose previously used value, select ... option.

"P" is the area of a single sample. It is possible to choose standard values ($0,02$ or $0,1\text{g/m}^2$) or enter specific value ("A" option).

To enter the value use the following keys:

- \blacktriangledown - digit increase,
- $\square \rightarrow$ - decimal point,
- $\rightarrow T \leftarrow$ - next digit,
- MENU* - end.

The result of paperweight measurement is finished with "=" mark pointing g/m^2 unit.

The balance is ready for the next measurements


Note:

1. Err-3 communicates that a sample was not put on the pan. The same communicate appears if single piece weight is less than one reading unit

15. Maintenance and small repairs

1. Keep the balance clean.
2. Take care, for during the operation some impurities may enter under balance base and into sensors vicinity. In case of impurities finding remove them.
3. In case of improper operation due to short power failure, switch off the balance by unplugging the power cable from the socket, and plug it again after a few seconds.
4. The message “Err-b”, displayed after switching on the unloaded balance, indicates the balance sensor mechanical damage.
5. All repairs must be performed by the authorized personnel.
6. To perform the repair contact the nearest service point.

Emergency messages:

Message	Reason	Recommendation
<i>C-1 ... 6 (over 1min.)</i>	autotest negative result	contact the service
<i>Err-b</i>	balance loaded during turning on	remove load from the balance
	balance sensor mechanical failure	contact the service
<i>l</i>	balance mechanical failure	contact the service
<i>H</i>	balance overload	remove load from the balance
	balance mechanical failure	contact the service
indicator does not work 	balance unstable, base vibrations, air blows	place the balance in location, assuring the indications stability
	balance damage	contact the service
-----	tare setting not finished	contact the service
- -	tare setting unsuccessful (too low load or B/G pressed)	zero the balance or pres B/G again
- -	zeroing with too high load	set the tare

Declaration of Conformity

We:

AXIS Spółka z o.o. 80-125 Gdańsk, ul.Kartuska 375B declare with full responsibility, that the
balances:

B60P

marked with CE mark are consistent with:

1. Standard PN-EN 61010-1:2004 Safety requirements for electrical equipment for measurement, control and laboratory use. General requirements harmonized with the directive 73/23/EEC (Low Voltage Directive),
2. Standard PN-EN 55022:2000 Electromagnetic compatibility (EMC) – Information technology equipment – Radio disturbance characteristics – Limits and methods of measurement, and PN-IEC 61000-4-3 Electromagnetic compatibility (EMC) – Part 4-3: Testing and measurement techniques - Radio-frequency electromagnetic field immunity test, harmonized with 89/336/EEC directive (electromagnetic compatibility).

Moreover, balances with the following information on rating plates:

- sticker with two digits of year of conformity assessment and with number of notified unit, performing the assessment,
- green metrological marking M,
- protective mark from notified unit



are made according to certificate of approval EC no. PL 04 020 and have EC verification, confirming conformity with:

3. Standard PN-EN 45501 Specification for metrological aspects of non-automatic weighing instruments, issued in December 1999, harmonized with the 90/384/EEC council directive, changed by 93/68/EEC council directive).

Additional information

- Conformity evaluation for the Council Directive 73/23/EEG and 89/336/EEG were carried out by Laboratorium Badawcze Oddziału Instytutu Elektrotechniki w Gdańsku, accredited to PCA.
- Approval certificate no. PL 04 020, issued by Central Office Of Measures in Warsaw (Notified Unit No. 1440).

Gdańsk, 9th of January, 2008

Authorized by the Director AXIS Sp. z o.o.:

Production manager mgr inż. Jan Kończak

A handwritten signature in black ink, appearing to read 'M. Kończak', is written over a horizontal line. Below the line, the word 'Signature' is printed.

NOTES